



The versatile SPOT™ Flex camera is a multipurpose workhorse camera with industry-leading performance specifications. Ultra-high resolution options up to 64 Mp are supported, along with full color sampling of the large format CCD sensor, providing accurate color rendition, crisp detail, and an incredible field of view. Cooling and 14-bit digital output reduce noise for low light level imaging. The suggested applications for this camera include widefield color imaging, high resolution documentation, fluorescence imaging, slide scanning, inspection, and phase contrast studies.

The SPOT Flex includes SPOT Basic software (Windows® and Mac® OS) which features real time image previewing, scale bar addition, report generation, time lapse, and measuring capabilities. IEEE 1394 Firewire™ connectivity enables effortless installation on computers and laptop portability. SPOT Flex cameras come with a 2-year warranty and are supported by an experienced team of imaging specialists.



# Flex™

4.0 Mp | COLOR

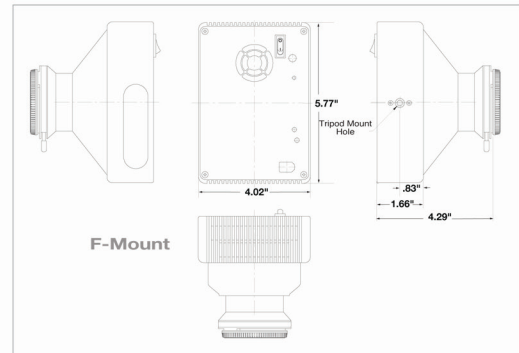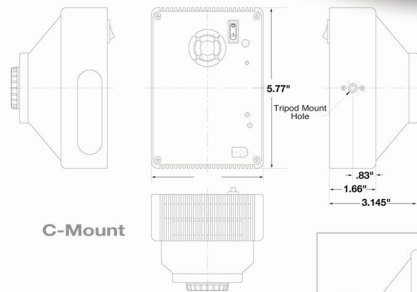

Features	Benefits
40 MHz live mode (dual channel 20 MHz)	High speed imaging for real time viewing
4 Mp Color Mosaic to 16 Mp* Cosite Color Sampling	Multiple resolution options provide the ability to capture a wide range of samples
-20° C differential cooling	Reduces thermal noise for low light level imaging
14-bit x 20 MHz capture	Extra bit depth is ideal for image enhancement
15.2 mm x 15.2 mm imaging area 21.4 mm FN	Large field of view to closely match the field of view in the microscope eyepiece
FireWire® interface	Allows for effortless installation, hot-swapping, and laptop connectivity for greater ease of use

\*Equivalent to 64 Mp Color Mosaic Resolution




**SPOT Imaging Solutions™**  
A division of Diagnostic Instruments, Inc.

Sensor	
Resolution	4.0 Mp to 16.0 Mp Cosite color sampled
Sensor Type	Kodak KAI-4020-CM with cover glass Color mosaic progressive scan interline CCD
Active Pixels	2048 x 2048
Pixel Size	7.4µm x 7.4µm
Optical Format	4/3" (4:3)
Active Imager Size	15.16 mm x 15.16 mm 21.43 mm diagonal
Exposure Time	1 ms to 536 seconds
Binning	1x1, 2x2 (Monochrome) 4x4 (Monochrome)
Cooling	-20°C differential from ambient via thermoelectric cooler with hermetically sealed, fan cooled heat sink
Read Noise	15 e <sup>-</sup> rms
Dark current	0.103 e/p/s mean value
Full Well Depth	33,333 e <sup>-</sup>
Non-linearity	<1%
Bit Depth	14-bit
Sensor Readout Frequency	40 MHz

Camera	
Lens Mount	C-mount (F-mount available)
Computer Interface	FireWire® IEEE 1394a interface
Trigger Connectors	Camera ready (TTL) External trigger input (TTL) Exposure sync output (TTL)
Tripod Mount	1/4"-20 UNC
Weight	C-Mount 2.4lbs. (1.09kg) F-Mount 2.4lbs. (1.09kg)
Size	(C-Mount) 5.77" x 4.02" x 3.15" (147mm x 102mm x 80mm)  (F-Mount) 5.77" x 4.02" x 4.29" (147mm x 102mm x 109mm)
Power Requirements	100-240 VAC, 47/63Hz
Operating Environment	0 to 30°C ambient, 0 to 80% relative humidity (non-condensing)
Certifications	CE, FCC Class A, EN60950 RoHS

### Quantum Efficiency

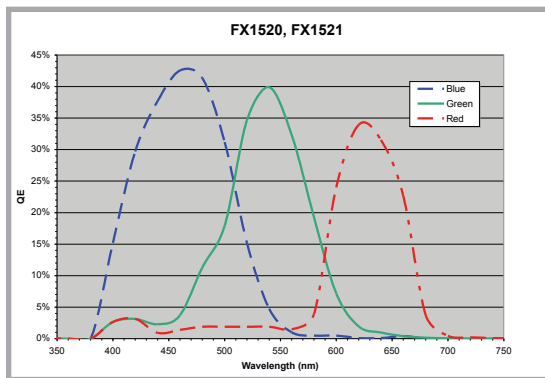

### Captured Frame Rates (Single Shot Mode)

Captured Frames per Second\*

Binning level	Region of interest			
	2048 x 2048	1024 x 1024	512 x 512	256 x 256
1 x 1	3.4	6.2	10.5	15.9
2 x 2 (Monochrome)	6.1	10.3	15.8	21.3
4 x 4 (Monochrome)	10.2	158.4	20.6	24.9

\*1 ms exposure, 20 Mhz readout with post-processing deferred, taken with 966 Mhz PIII processor running Windows XP. Single shot mode. Capture rates on other computers OS platforms may vary.

SPOT Software	
SPOT Basic Software Features	Color Live mode viewing window with controls Image capture window Image capture setups are user defined, named and saved Auto or manual exposure for Live and captured images Zoom control on Live and captured images Auto white balance Spot or average metering Flatfield correction Annotation Calibration mark Measurement tools Image enhancement tools Sequential image capture and playback Report generator Interactive print dialog Help menu 3rd Party Software Interface: Standard Twain for Windows®, SPOT Image Capture Application & AppleEvent for Mac®
Computer Operating System Compatibility	Windows®: XP SP3 (32-bit), Vista (32 or 64-bit), 7 (32 or 64-bit) and 8 (32 or 64-bit) Mac®: OS 10.6 - 10.8
Minimum Computer Requirements	PC: 1GHz Pentium 4 (or equivalent x86 SSE2 processor), 1 GB RAM, IEEE 1394a port Mac®: Intel processor, 1 GB RAM, IEEE 1394a port Mac® and PC video: 24-bit color @ 1024 x 768 resolution
SPOT Advanced Software	Available (includes macro feature, floating taskbars, custom dialogs and image archiving database)
SPOT Advanced Software Modules (Optional)	Available (adds image display intensity scaling, peripheral device control, extended depth of focus and extended third party interface)
Software Development Kit	Available for Mac® and Windows®
Native Drivers for 3rd Party Software	Available (call or visit our website)

Items Included	
Camera Head	With lens cap
FireWire® cables	Data cables (6-pin and 6 to 9 pin)
Power supply	With power cord, 2.29" (58mm) x 5.20" (132mm) x 1.18" (30mm), 0.70 lbs. (0.32kg)
Quick Start Installation Guide	Instructions on how to get up and running quickly
Software CD	Includes SPOT Basic software and User Guide
Warranty	2 years from date of purchase

Catalog Number	C-Mount	F-Mount
Flex Color	FX1520	FX1521



Specifications are typical and subject to change.  
**SPOT Imaging Solutions™**  
 A division of Diagnostic Instruments, Inc.